


X-RAY FILM DEVELOPING TANK

USER MANUAL

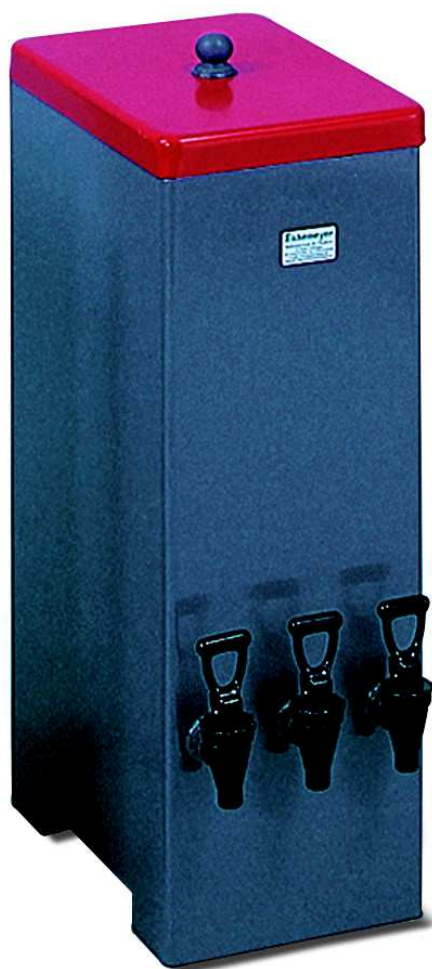

Art. Nr. 706920
Art. Nr. 706922

TELEPHONE +49 7461 96 580 0

www.eickemeyer.com

veterinary technology for life
Eickemeyer®

Contents

Contents 2

1 Components of Developing Tank..... 3

2 Product Function and Operation 3

3 Cleaning and maintenance..... 3

1 Components of Developing Tank

The developer tank has 3 compartments:

The middle for water rinse.

The left and right for the developer and fixer.

2 Product Function and Operation

Prior to the first use please clean/rinse all compartments with warm water.

Check valves for leaks by filling again lukewarm water up to $\frac{3}{4}$ of the maximum filling level.

If no leaks occur, you can now fill in both 1.8 litre concentrate in the two external tanks (applies to concentrate which must be diluted 1:4)

Stir the contents thoroughly.

If you have a developer tank with heater, you can adjust the temperature at the back of the heating element.

Set the setting pin about the centre position, which is about 26 ° C water temperature.

Note: Due to the small internal dimensions of the tank (9 l volume) always only one film should be processed

- Recommended time in the developer at 25/26° C: 2 min
- Please fill in the amount of loss by spreading of the developer tank with developers concentrate.

.

3 Cleaning and maintenance

After about 4 weeks or 20 films, the content should be replaced. We recommend to exchange the water after about 3-4 films.

Before each new batch of fresh concentrate you have to observe that the inside of the developer tanks should be thoroughly cleaned. Rinse with clean water.

GERMANY

EICKEMEYER KG
Eltastraße 8
78532 Tuttlingen
T +49 7461 96 580 0
F +49 7461 96 580 90
E info@eickemeyer.de
www.eickemeyer.de

ITALY

EICKEMEYER S.R.L.
Via G. Verdi, 8
65015 Montesilvano (PE)
T +39 085 935 4078
F +39 085 935 9471
E info@eickemeyer.it
www.eickemeyer.it

CANADA

EICKEMEYER Inc.
250 Briarhill Dr.
Stratford, Ont. Canada
N5A 7S2
T +1 519 273 5558
F +1 519 271 7114
E info@eickemeyervet.ca
www.eickemeyervet.ca

NETHERLANDS

EICKEMEYER B.V.
Bedrijventerrein
Pavijen-West Bellweg 44
4104 BJ Culemborg
T +31 345 58 9400
F +31 345 58 9401
E info@eickemeyer.nl
www.eickemeyer.nl

SWITZERLAND

EICKEMEYER AG
Sandgrube 29
9050 Appenzell
T +41 71 788 23 13
F +41 71 788 23 14
E info@eickemeyer.ch
www.eickemeyer.ch

DENMARK

EICKEMEYER ApS
Lysbjergvej 6, Hammelev
6500 Vojens
T +45 7020 5019
F +45 7353 5019
E info@eickemeyer.dk
www.eickemeyer.dk

POLAND

EICKEMEYER Sp. z o.o.
Al. Jana Pawła II 27
00-867 Warszawa
T +48 22 185 55 76
F +48 22 185 59 40
E info@eickemeyer.pl
www.eickemeyer.pl

UNITED KINGDOM

EICKEMEYER Ltd.
3 Windmill Business Village
Brooklands Close
Sunbury-on-Thames
Surrey, TW16 7DY
T +44 20 8891 2007
F +44 20 8891 2686
E info@eickemeyer.co.uk
www.eickemeyer.co.uk